




IRD120 Infrared Thermography and Disinfection Channel DATA SHEET

PRODUCT DETAILS

Product name:	Intelligent infrared thermography and disinfection channel with fever alarm system		
Working theory:	<p>1, It uses infrared thermal imaging non-contact temperature measurement, auto focus, 1 second temperature measurement, thermal imaging accuracy can reach $\pm 0.3^{\circ}\text{C}$, to avoid cross infection caused by contact temperature measurement.</p> <p>2, Thermal disinfectant spray nozzle with infrared sensor system in the channel sterilizes the body automatically surface at 360 degrees without dead angle, and only needs to stay for 5 seconds to kill bacteria and viruses.</p>		
Features:	<ul style="list-style-type: none"> • Non-contact: Automatic temperature measurement with thermography tech, which reduces the risk of cross infection. • Auto-focus, temp measurement in one second: precisely up to $\pm 0.3^{\circ}\text{C}$. • Auto fever alarm: Auto alarm with LED flash and voice broadcast for people with abnormal temperature, which reduces the workloads of the operator. • 360° disinfection: Induction spray of disinfection in 5 seconds. • Voice broadcast: Auto voice broadcast of detected person temperature. • Software for report analysis, the information is traceable: Do not miss any suspected person, each person picture and thermal image and temperature data will be recorded automatically for convenient tracking. 		
Item No.:	IRD120	Type:	Vertical and movable
Technology:	Auto thermography and disinfection with infrared sensor	Power supply:	Plug-in
Certificate:	ROHS, will get CE 10 days later	Screen:	Touch screen
Functions:	Temperature measurement, disinfection, fever alarm, data record and download	Screen size:	12.1 inches
Optional functions:	Hand induction disinfectant and alcohol concentration detection	Equipment size:	100x100x200cm
Detection Accuracy:	$\pm 0.3^{\circ}\text{C}/\pm 0.5^{\circ}\text{C}$	System:	Android 7.0
working temperature:	0-45 $^{\circ}\text{C}$.	Detection distance:	1.5m
Capacity of photos:	50,000pcs	Lack of liquid alarm:	Yes
Fogging quantity:	20L (1, consume 8ml each time 2, can be used for 2500 people)	Scope of application:	Community/hospital/school /mall/hotel/cinema/enterprise/factory/park
Port of loading:	Wuhan, China or any port you required	Personal information (Name, and phone) record:	No
Machine Weight:	100KG	Input voltage:	100-240V , 50/60HZ
Input power:	200W	Plug type:	All the plug type is available.

PHOTO DETAILS

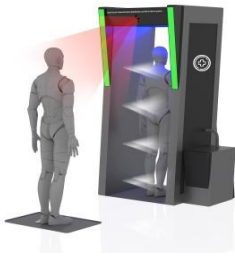

PACKAGING INFO

N.W. of machine:	100KG	
Machine size:	100x100x200cm	
Sales packaging:	3 boxes	(70*70*40cm, 20*110*250cm, 100*100*40cm)
ETD: 12 working days	G.W.: 180KG	
ETA: About 7-10 working days after the shipment by Air.		
Tracking No.:	Available, will be advised after the shipment.	
Payment item:	T/T etc.	
Warranty period:	One year	